



TISSUE EXTRACTED PROTEIN PROTECTION FROM ENDOGENOUS PROTEASES USING PROTEIN SAFE PRECELLYS LYING KIT

/ CONTEXT

Homogenization of biological samples for protein extraction using Precellys® Evolution and Cryolys® Evolution ensures a high yield of protein recovery and prevents protein from heat degradation during this process. However, all proteins, including endogenous proteases, being protected during this step, a subsequent degradation of the proteins of interest can occur during and after the homogenization. Therefore, it is highly advised to protect proteins of interest from degradation by endogenous proteases released by the sample.

In this Application Note, we highlight how the use of the Protein Safe Precellys® Lysing kits ensures an efficient protection of samples from endogenous proteases in different kinds of tissue homogenates.

/ MATERIALS

- Precellys® Evolution and Cryolys® Evolution
- Precellys Lysing kits CK14 (ref. P000912-LYSK0-A.0) and Protein Safe Precellys® Lysing kit CK14 (ref. P000973-LYSK0-A.0)
- Rat organs: liver, brain, heart, lung from Janvier labs
- Assay buffer (Potassium phosphate base buffer)
- Fluorescent substrate of proteases: universal protease substrate resorufin labelled (Sigma 11734334001)
- Fluorimeter PerkinElmer Type LS50B

/ PROTOCOL

The different tissue types (liver, brain, heart and lung) were cut off in-pieces. 50mg of sample were then loaded in either CK14 Protein Safe lysing tubes or in regular CK14 Lysing tubes. Volume was then completed with 1,6mL of buffer. Tubes were placed on the Precellys® Evolution equipped with Cryolys® Evolution instrument and processed using a generic homogenization program with the following settings:

- 2 x 30sec at 6500rpm
- Break 15sec
- Cooling: 4°C.





400µl of the obtained homogenate were then transferred in a new tube to proceed to proteases activity assay. Assay is based on incubation of the homogenate with a universal proteases resorufin labelled substrate. Active proteases degrade this substrate and release a fluorescent molecule. The level of fluorescence detected on the sample by fluorimeter is then related to protease activity. In case proteases are inactive (Protein Safe inhibitor), the reaction cannot occur and no or low fluorescence is observed.



TISSUE EXTRACTED PROTEIN PROTECTION FROM ENDOGENOUS PROTEASES USING PROTEIN SAFE PRECELLYS LYSING KIT

/ RESULTS

Fluorescence levels reflecting the proteases activity levels were measured in the different samples following the method described previously. The different activity levels measured for a same sample extracted with Classical Precellys® Lysing kit or the Protein Safe Lysing kit are represented in the figures 1 to 4 for Heart, Liver, Lung and Brain.

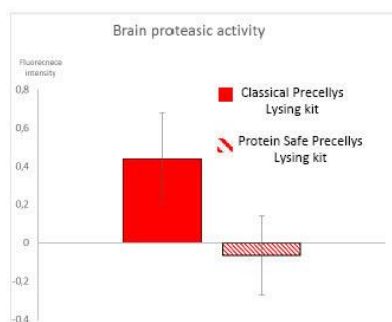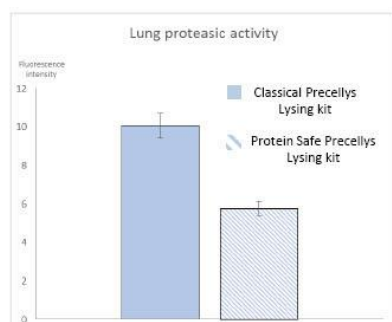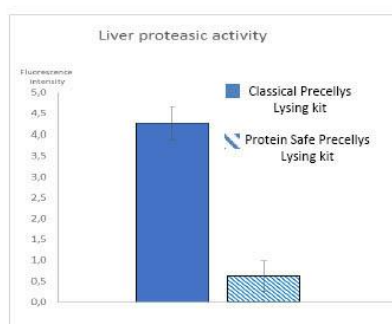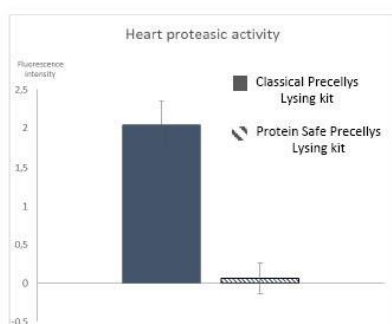

The protease activity significantly decreases when samples are homogenized using the Protein Safe Precellys® Lysing kits.

In the case of Heart, Liver and Brain, the inhibition of proteases is close to 100% and in the case of Lung, close to 50%. The difference of efficiency can be explained by the expression levels of proteases in the different organs or by the nature of the proteases themselves, since Metalloproteases are not inhibited by the Protein Safe buffer (for MMP EDTA can be added to the tube but has not been evaluated in the present case).

/ CUSTOMER



PREC-026-DU140-A

/ CONCLUSION

The Protein Safe Precellys® lysing kit is efficient to inhibit proteases activity in various organs at a significant level during and after the homogenization process with Precellys® and Cryolys® Evolution. Proteins of interest are then protected from endogenous degradation once the homogenization step finished.



□東京 〒162-0805 東京都新宿区矢来町 113 番地 TEL (03)3235-0661(代) / FAX (03)3235-0669

□大阪 〒532-0005 大阪市淀川区三国本町2丁目12番4号 TEL (06)6396-0501(代) / FAX (06)6396-0508

※会社名および商品名は、各会社の商標または登録商標です。
※本カタログに記載の規格・仕様・外観は予告なく変更する場合がありますので御了承下さい。

<http://www.technosaurus.co.jp>

【Precellys Evolution でのタンパク抽出において、PROTEIN SAFE LYSING KIT を用いることによる
試料内在性タンパク分解酵素抑制効果の検証】

Precellys-APS-1902-106